

ULTRASENSITIVE THYROID STIMULATING HORMONE (U-TSH) ENZYME IMMUNOASSAY TEST KIT

Catalog Number: BC-1003



mymsupply.com

Enzyme Immunoassay for the Ultrasensitive Quantitative Determination of Thyroid Stimulating Hormone (U-TSH) in Human Serum

FOR IN VITRO DIAGNOSTIC USE ONLY

Store at 2 to 8°C.

PROPRIETARY AND COMMON NAMES

Ultrasensitive TSH Enzyme Immunoassay

INTENDED USE

For the quantitative determination of the thyroid stimulating hormone (TSH) concentration in human serum.

INTRODUCTION

The determination of serum or plasma levels of thyroid stimulating hormone (TSH or thyrotropin) is recognized as a sensitive method in the diagnosis of primary and secondary hypothyroidism.¹ TSH is secreted by the anterior lobe of the pituitary gland and induces the production and release of thyroxine and triiodothyronine from the thyroid gland.² It is a glycoprotein with a molecular weight of approximately 28,000 daltons, consisting of two chemically different subunits, alpha and beta.³

Although the concentration of TSH in the blood is extremely low, it is essential for the maintenance of normal thyroid function. The release of TSH is regulated by a TSH-releasing hormone (TRH) produced by the hypothalamus. The levels of TSH and TRH are inversely related to the level of thyroid hormone. When there is a high level of thyroid hormone in the blood, less TRH is released by the hypothalamus, so less TSH is secreted by the pituitary. The opposite action will occur when there is decreased thyroid hormone in the blood. This process is known as a negative feedback mechanism and is responsible for maintaining the proper blood levels of these hormones.^{4,5}

TSH and the pituitary glycoproteins: luteinizing hormone (LH), follicle-stimulating hormone (FSH), and human chorionic gonadotropin (hCG), have identical alpha chains. The beta chains are distinct but do contain regions with identical amino acid sequences. These regions of homology can cause considerable cross-reactivity with some polyclonal TSH antisera. The use of a monoclonal antibody in this TSH ELISA test eliminates such cross-reactivity, which could result in falsely elevated TSH values in either menopausal or pregnant females -- a population whose evaluation of thyroid status is clinically significant.^{6,7,8}

PRINCIPLE OF THE TEST

The Ultrasensitive TSH ELISA test is based on the principle of a solid phase enzyme-linked immunosorbent assay.^{9,10} The assay system utilizes a unique monoclonal antibody directed against a distinct antigenic determinant on the intact TSH molecule. Mouse monoclonal anti-TSH antibody is used for solid phase immobilization (on the microtiter wells). A goat anti-TSH antibody is in the antibody-enzyme (horseradish peroxidase) conjugate solution. The test

sample is allowed to react simultaneously with the two antibodies, resulting in the TSH molecules being sandwiched between the solid phase and enzyme-linked antibodies. After a 2-hour incubation at room temperature, the wells are washed with water to remove unbound labeled antibodies. A solution of TMB Reagent is added and incubated for 20 minutes, resulting in the development of a blue color. The color development is stopped with the addition of Stop Solution changing the color to yellow. The concentration of TSH is directly proportional to the color intensity of the test sample. Absorbance is measured spectrophotometrically at 450 nm.

REAGENTS

Materials provided with the kit:

- Murine Monoclonal Anti-TSH-coated microtiter wells.
- Set of Reference Standards: 0, 0.1, 0.5, 2, 5, and 10 μ IU/ml, lyophilized.
- Enzyme Conjugate Reagent, 13 ml.
- TMB Reagent (One-Step), 11 ml.
- Stop Solution (1N HCl), 11 ml.

Materials required but not provided:

- Equipment capable of shaking microplates at 175 RPM.
- Precision pipettes: 100 μ l and 1.0 ml.
- Disposable pipette tips.
- Distilled water.
- Vortex mixer or equivalent.
- Absorbent paper or paper towel.
- Graph paper.
- Microtiter plate reader.

SPECIMEN COLLECTION AND PREPARATION

Serum should be prepared from a whole blood specimen obtained by acceptable medical techniques. This kit is for use with serum samples without additives only.

STORAGE OF TEST KIT AND INSTRUMENTATION

Unopened test kits should be stored at 2-8°C upon receipt and the microtiter plate should be kept in a sealed bag with desiccants to minimize exposure to damp air. Opened test kits will remain stable until the expiration date shown, provided it is stored as described above. A microtiter plate reader with a bandwidth of 10nm or less and an optical density range of 0-2 OD or greater at 450 nm wavelength is acceptable for use in absorbance measurement.

REAGENT PREPARATION

1. All reagents should be brought to room temperature (18-25°C) before use.
2. Reconstitute each lyophilized standard with 1.0 ml distilled water. Allow the reconstituted material to stand for at least 20 minutes and mix gently. Reconstituted standards will be stable for up to 30 days when stored sealed at 2-8°C

0	0.009
0.1	0.050
0.5	0.218
2	0.801
5	1.800
10	3.191

ASSAY PROCEDURE


1. Secure the desired number of coated wells in the holder.
2. Dispense 100 μ l of standards, specimens, and controls into appropriate wells.
3. Dispense 100 μ l of Enzyme Conjugate Reagent into each well.
4. Thoroughly mix for 30 seconds. It is very important to mix them completely.
5. Incubate at room temperature (18-25°C) with *shaking at 175 RPM, for 120 minutes*.
6. Remove the incubation mixture by flicking plate contents into a waste container.
7. Rinse and flick the microtiter wells 5 times with distilled or deionized water. (*Please do not use tap water.*)
8. Strike the wells sharply onto absorbent paper or paper towels to remove all residual water droplets.
9. Dispense 100 μ l of TMB Reagent into each well. Gently mix for 10 seconds.
10. Incubate at room temperature for 20 minutes.
11. Stop the reaction by adding 100 μ l of Stop Solution to each well.
12. Gently mix for 30 seconds. *It is important to make sure that all the blue color changes to yellow color completely.*
13. Read absorbance at 450 nm with a microtiter well reader within 15 minutes.

CALCULATION OF RESULTS

1. Calculate the mean absorbance value (A_{450}) for each set of reference standards, controls and patient samples.
2. Construct a standard curve by plotting the mean absorbance obtained from each reference standard against its concentration in μ U/ml on graph paper, with absorbance values on the vertical or Y axis, and concentrations on the horizontal or X axis.
3. Use the mean absorbance values for each specimen to determine the corresponding concentration of TSH in μ U/ml from the standard curve.

EXAMPLE OF STANDARD CURVE

Results of a typical standard run with absorbency readings at 450 nm shown in the Y axis against TSH concentrations shown in the X axis. This standard curve is for the purpose of illustration only, and should not be used to calculate unknowns. Each user should obtain his or her own data and standard curve.

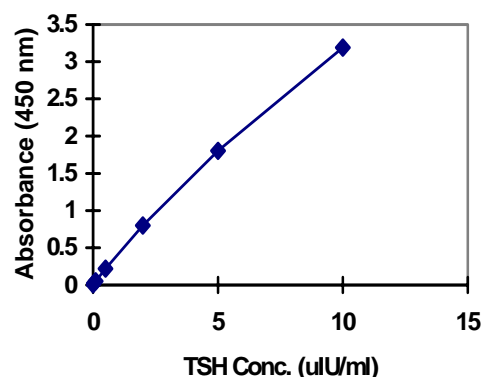

EXPECTED VALUES AND SENSITIVITY

The mean TSH values based on 160 random normal adult blood samples, is 1.6 (0.4-7.0) μ U/ml. TSH levels exceeding 10 μ U/ml, suggest primary hypothyroidism. Low or undetectable TSH levels may be normal, but may also indicate secondary hypothyroidism (insufficient secretion of TSH or TRH). Low levels may also be due to hyper-secretion of T-3 and T-4 due to Grave's disease or thyroiditis. Differential diagnosis is best achieved by simultaneous determination of TSH and free T-4 levels in serum.

The minimum detectable concentration of TSH by this assay is estimated to be 0.05 μ U/ml.

LIMITATIONS OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the package insert instructions and with adherence to good laboratory practice.
2. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.
3. Serum samples demonstrating gross lipemia, gross hemolysis, or turbidity should not be used with this test.
4. The results obtained from the use of this kit should be used only as an adjunct to other diagnostic procedures and information available to the physician.

REFERENCES

- ¹ Burger, H. G., Patel, Y. C., Thyrotropin releasing hormone -TSH Clinic. Endocrinol. and Metab., 6, 831-00(1977).
- ² Ezrin, C., The Thyroid, S. C. Werner and S. H. Ingbar (eds.), Harper and Row, Hagerstown, MD, 9, 174-178 (1978).
- ³ Pierce, J. G., Endocrinology, 89, 1331-1344 (1971).
- ⁴ Berger, S. and Quinn, J. L., Fund. Clin. Chem., N. W. Tietz (ed.), W. B. Saunders Co., Phila., PA 14, 824-848 (1976).
- ⁵ Utiger, R. D., The Thyroid, S.C. Werner and S. H. Ingbar (eds.), Harper and Row, Hagerstown, MD, 9, 196-205 (1978).

TSH (μ U/ml)	Absorbance (450 nm)
-------------------	---------------------

- ⁶ Soos, M. and Siddle, K., J. Immun. Methods, 51, 57-68 (1982).
- ⁷ Wada, H. G., Danisch, R. J., Baxter, S. R., et al, Clin. Chem.,28, 1862-1866 (1982).
- ⁸ Snyder, P. J. and Utiger, R. D., J. Clin. Endocrinol. Metab., 34, (1972).
- ⁹ Engall, E., Methods in Enzymology, Volume 70, Van Vunakis, H. and Langone, J. J. (eds.), Academic Press, New York, 419-492(1980).

040203



DISTRIBUOTOR:
MYM Laboratory & Medical Supply, Inc.
8684 Ave. de la Fuente Ste. 14
San Diego CA. 92154
Tel. (619) 710-0126 Fax. (619) 710-0297
www.mymSupply.com
email: mail@mymSupply.com